



BEAUCARE[®]

Better Value. Personal Service. Trusted Partner.


ProBed
Range

ProBed Byland Bariatric Profiling Bed Instruction Manual

PRO-BEDX4



USER MANUAL

www.beucare.com | sales@beucare.com | 01423 873666

1.0 Technical Specifications

1.1 Safe working load

The maximum safe working load (SWL) for this bed is 200kg / 31.5 stones

1.2 Weight of the Bed.

The weight of the bed is 120 kg. This weight can vary according on the choices of accessories installed. For easy handling, transport and storage the bed is divided into four parts. See point 2 for instructions on how to assemble and dismantle the bed.

1.3 Accessories

Features:

- Hand controller
- Castors with individual braking system
- Wooden headboard
- Wooden footboard
- Lower bedframe sub assembly
- Upper bedframe sub assembly
- Wooden side rails

1.4 Storage.

The bed should be kept in a dry and ventilated place.

1.5 Bed relocation

Before you move the bed, it must be on its lowest position. Always use the headboard or the footboard to move it. Never use the side rails to move the bed. Always handle the bed carefully. Apply the brakes when the bed is in place. Remembering to remove the mains power supply and ensuring all cables are positioned so they do not interfere with the castors or persons moving the bed.

1.6 Measurements.

Mattress measurements suitable for this bed: 1100 x 1950 mm.



Maximum angle: Backrest: 63°.

Footboard: 30°.

Trendelenburg (headboard section descent) 12°.

Against trendelenburg (foot section descent) 12°.

Total length of bed: 2000 mm.

Overall width: 1200 mm dimensions of pedestal: 90 cm x 1950 cm.


Bed height without mattress (with 75mm diameter castors).

Highest position: 800 mm-lowest position: 200 mm

Note: the dimensions specified in this manual do not consider the manufacturing tolerances.

2.0 Bed Installation

The assembly process can be done by one or two people. However due to the size of the structures, we recommend two people for the assembly to be more secure.

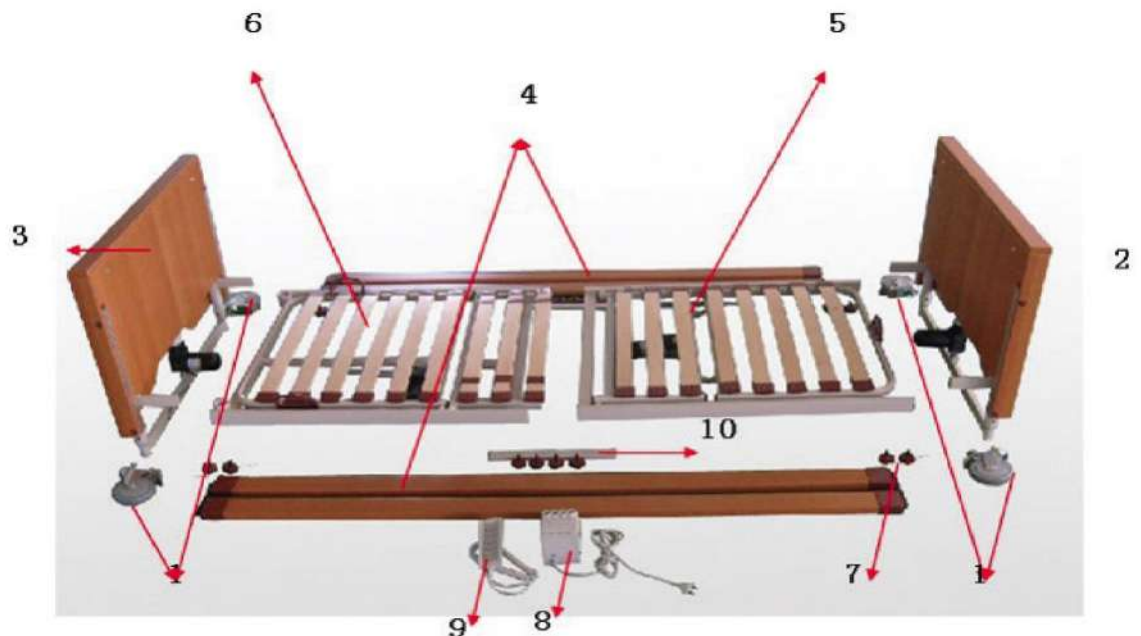

Item 1	Name	Description	Qty
1	Castors	75mm wheel, with individual brakes	4
2	Headboard and 'H' bed lifter subassembly	Assembled headboard and 'H' bed lifter support and base with lifted motor	1
3	Footboard and 'H' bed lifter subassembly	Assembled footboard and 'H' bed lifter support and base with lifted motor	1
4	Side rails	Four wood side rails with plastic cover and springs	4
5	Backrest bed frame subassembly	Assembled bed frame / backrest, containing back – rest motor	1
6	Lower knee-break bed frame subassembly	Assembled bed frame / knee-break, with knee-break motor	1
7	Fixing screws	With plastic grips	4
8	Linak control box		1
9	Linak hand controller		1
10	Steel Connecting Bars	Bars to join 5 & 6 subassembly's	2

2.2 Installation process

(A) Install Head/Footboard (2&3) to bed frame subassembly (5&6): Stand Install Head/Footboard (2&3) up against a wall with the brakes applied. (2) and (3) are the same so either can be used either end. Open the spring locking mechanism to allow bed frame to be applied as shown.



Then located the metal lugs of the bed frame into the receiving notches on the bed end and push down firmly so they sink to the bottom.



Take care to close the spring locking mechanism into its locked position.



(B) Connect Bed frame Backrest subassembly (5) to (6) Hold up the foot end bed frame subassembly, and put the steel connection bars (10) into tube of bed frame. Use the fixing screws (7) to secure in, hand tighten.


The steel connection bars are located on the transport frame on delivery.



Please note: tighten the screws.

(C) Connect Bed frame Foot end subassembly (6) to (5) The assembly method is same as 2.1.2 (c) described above.

2.3 Connect plus for motors

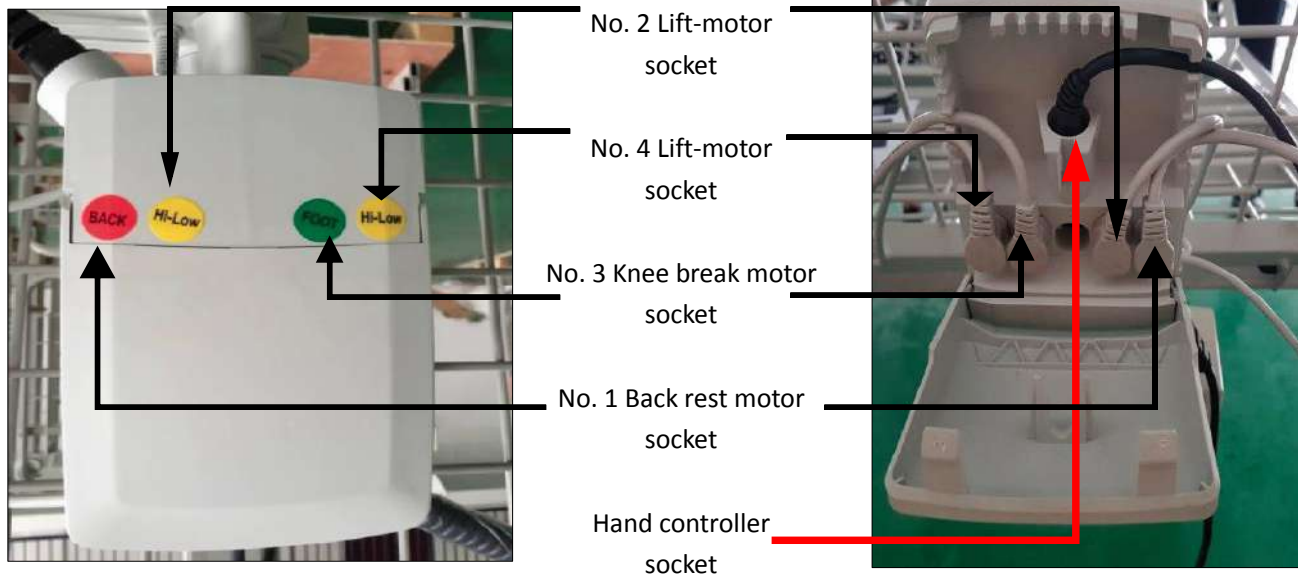

The control box (8) is already installed on the bed frame subassembly (6). Remove the cap, locate the plugs for the 4 motors and the hand controller (9) and insert them firmly into the correct sockets as shown in the above image.



Ensure the cables are routed through the metal holders located on the bed Platform




Ensure the safety key is used to turn on the required controls on the handset.

2.4 Side rails (4) assembly

Next plug the power cable into the mains and switch on at the wall. Raise bed up to allow side rails to be installed.

Side Rail Finger Slider Components:

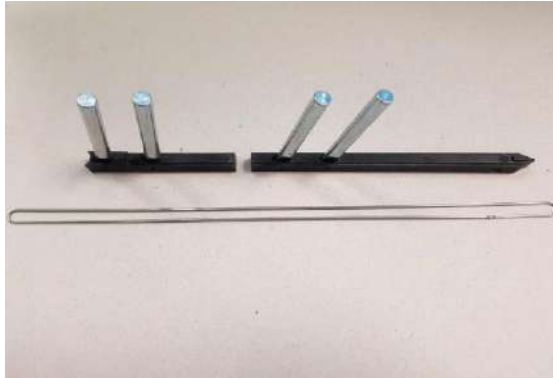

- A) Assemble finder slider components as in Fig.1.
- B) Slide the two plastic sections together as shown in Fig.2.
- C) Ensure wire loop is installed centrally in the curved groove provided Fig.2 and 3.
- D) Depress the lower release button on one end and slide the unit into the channel past the lower release button so it is held in the channel as Fig.4.
- E) Repeat A and B above with the 2nd slider and insert into the two side rails as shown in Fig.5 and 6.
- F) Fix the opposite end of the side rails into the fingers inserted into the channel in part C above
- G) Now complete the installation by sliding the 2nd slider into the channel by depressing the lower release button as in C above.
- H) Pull the rails up until they click past the top release buttons on both sides.

How the sliders work in the channel:



Fig.1

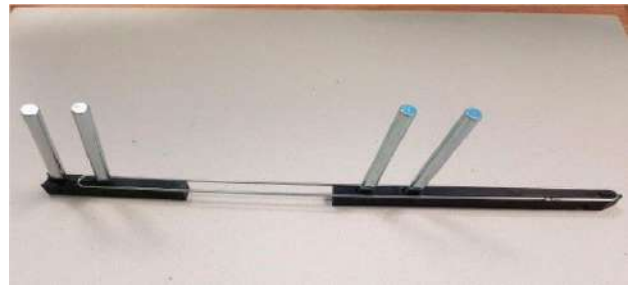


Fig. 2



Fig.3



Fig.4



Fig.5



Fig.6

3.0 After-sales service

If there are quality problems not caused by the user within a week after purchase, or quality problems, under normal use and care, our company ensures the free repair of the equipment once it has been purchased within one year . If there are problems after one year after the purchase, users should contact our after-sales department or store dealers with invoice to be sold the accessory or part.

The following cases are not included in free maintenance:

Damage that occur under the use.

Damage that occurs due to mounting and dismounting is made by users.

Damage caused by accidents or by the user.

4.0 Warranty

This warranty covers only those defects of parts and components of equipment described in normal use – according to the instructions of the operating manual accompanying the same, and which are provided by manufacturers.

This warranty will be automatically cancelled if the equipment undergoes repairs by unauthorised persons, receive ill-treatment or suffer damage from accidents, falls, voltage variations and overloading above the specified, or any unpredictable occurrence, arising from improper use of the equipment by the user.

The structure of this bed is covered by warranty for a period of 24 months. Where accessories and electrical part are covered for a period of 12 months warranty from the date of purchase.

Repairs covered by warranty must be carried out by our company, and the costs of transport owner's liability. The value of the repair.

5.0 Declaration of conformity

DECLARES: that the bed "decubitus electric" complies with the essential requirements laid down in annex I to Directive 93/42/EEC of 14 June, on your current writing and of Decree-Law No. 145/2009 of June, applicable, and does not commit the State medical or patient safety, nor even the safety and health of users or third parties, where appropriate, when used under the conditions and for the purposes intended.

Note: In case of repair should present this manual.

We hereby declare that the above mentioned product meets the provision of the Medical Device Regulation (EU) MDR 2017/745 of the European Parliament and of the Council of 5 April 2017 on medical devices. The EU declaration of conformity is issued under the sole responsibility of the manufacturer.



Standard Applied:

EN ISO 13485:2016	EN ISO 14971:2012	EN 60601-1:2006/A1:2013	EN 60601-1-2:2015
EN 60601-2-52:2010	EN 1041:2008	BS EN ISO 15223:1-2016	

Beaucare Medical Ltd
Hornbeam Square North
Hornbeam Park
Harrogate
North Yorkshire
HG2 8PA
Tel.: 01423 873666
www.beaucare.com
sales@beaucare.com

